


Tutorials



The Cleft Care app is designed to automatically identify hypernasality and speech acceptability in children with cleft lip and palate, and to refer those who need intervention to AIISH cleft care services for therapy or surgery. This app allows a Community Health Worker to record a speech sample during a patient visit. The speech sample will be analyzed by our automated system, and the results will be sent via email to the All India Institute of Speech and Hearing (AIISH), where speech pathologists will review the results and the samples uploaded to a secure server by the app and make recommendations for the child. If follow-up by the cleft palate team is recommended, the AIISH team will contact the CHW and the parent.

Instructions on how to use the application:

Login Screen: The CHW logs into the app using a secure login screen that verifies each login against a database of authorized CHW email addresses maintained by the AIISH speech team; any unauthorized email triggers an error message to prevent misuse. The app supports both English and Kannada, and users select their preferred language from a dropdown list so that all stimuli and interface text appear in the chosen language.

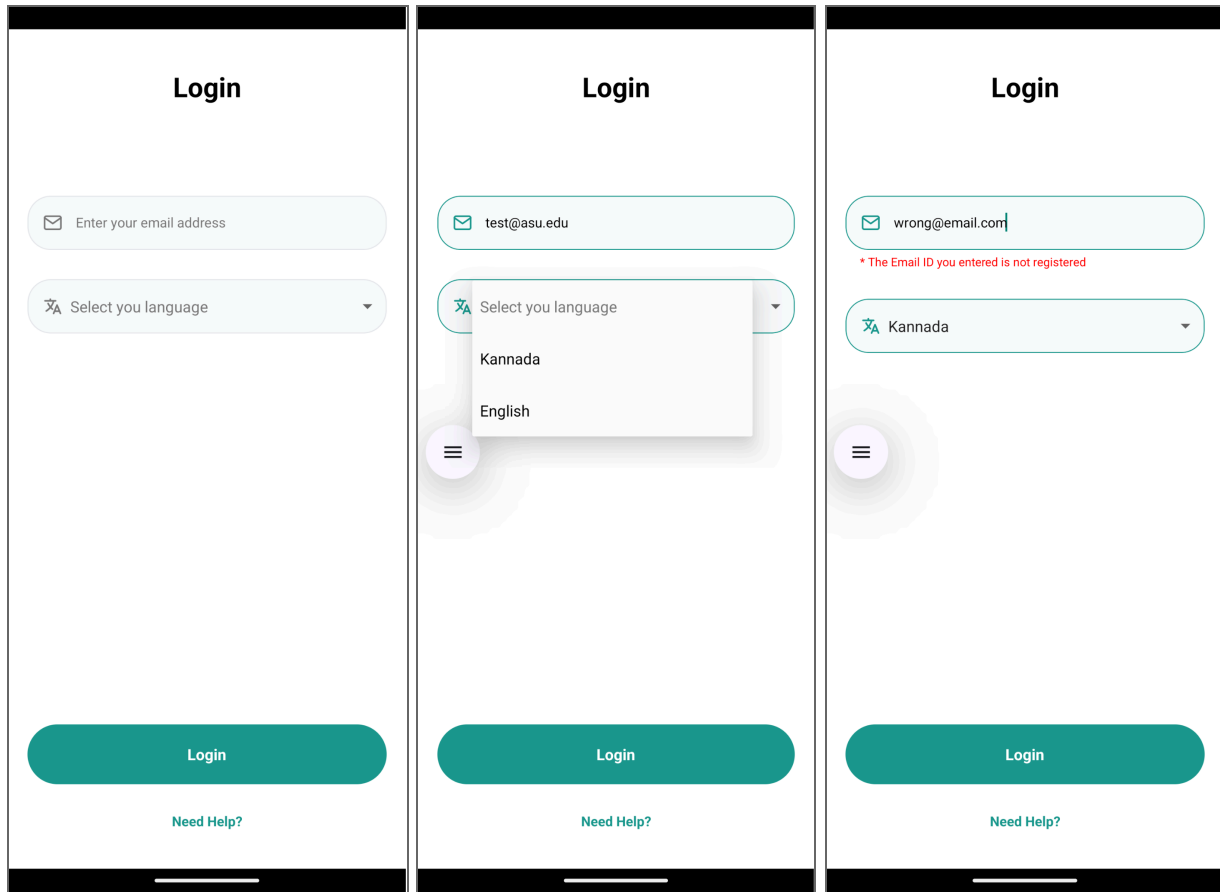

After login, three tabs appear at the bottom of the screen—Home, Tutorials, and the CHW Profile page—while the top-right corner includes three icons: a search tool to find any patient record by name or record ID, an “Add new record” button to register a new patient, and a help icon to send an email to the support center.

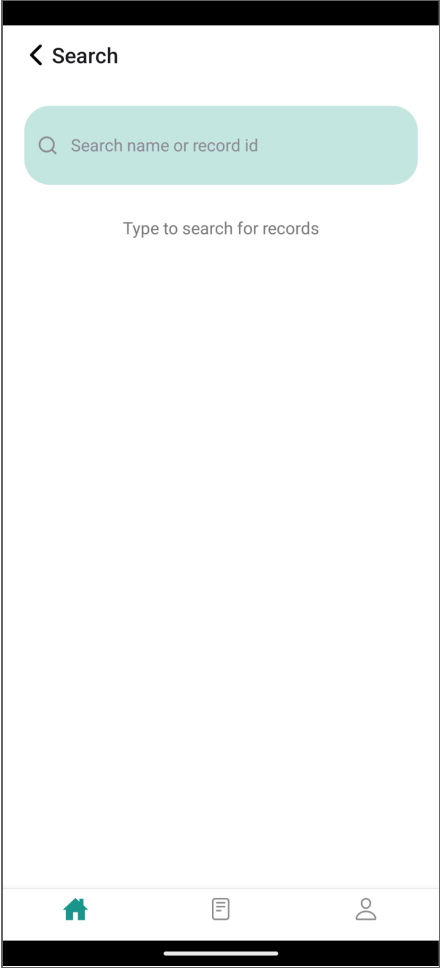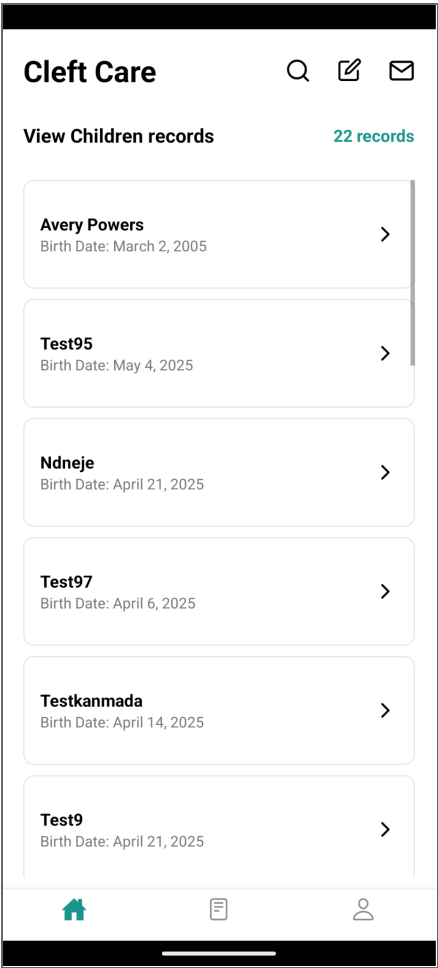

The image displays two side-by-side mobile application screens. The left screen, titled 'Add Children Record', features a back arrow and a title bar. Below the title, it instructs the user to 'Please enter child's information to continue with the audio recording'. The form includes several input fields: 'Enter your name', 'Select Birth Date', 'Select Gender' (with a dropdown arrow), 'Hearing Status' (with a dropdown arrow), 'Enter your Address', 'Contact number', and 'Take/Choose Photo (Optional)'. A 'Next' button is positioned at the bottom of the form. The right screen, titled 'Help Center', has a close button (X) in the top right corner. It prompts the user with 'Need help? Please let us know'. Below this is an email input field containing 'cleftcareasu@gmail.com' and a large text area for the message, labeled 'Enter your message here'. A 'Submit' button is located at the bottom of the screen. Both screens have a bottom navigation bar with three icons: a home icon, a document icon, and a user profile icon.


Adding a New Patient: To add a new patient, the CHW taps the “Add new record” button/icon and enters demographic details such as name, date of birth, gender, hearing status (presence or absence of hearing loss), residential address, and contact number, and can also upload a photo of the patient’s cleft palate as a medical record. After completing these fields, the CHW taps “Next” to review and complete the consent form, where they type the caregiver’s name and provide a signature to accept the terms. Once the consent form is signed, the CHW proceeds to the recording and begins recording the patient’s speech samples.

The image displays three sequential mobile app screens for a healthcare application.


- Screen 1: Add Children Record**
 - Header: < Add Children Record
 - Instruction: Please enter child's information to continue with the audio recording
 - Form fields:
 - Enter your name (with person icon)
 - Select Birth Date (with calendar icon)
 - Select Gender (with person icon and dropdown arrow)
 - Hearing Status (with smiley face icon and dropdown arrow)
 - Enter your Address (with location pin icon)
 - Contact number (with phone icon)
 - Take/Choose Photo (Optional) (with camera icon)
 - Next button (disabled)
- Screen 2: Consent Form**
 - Header: < Consent Form
 - Instruction: Please carefully read and review the consent form to ensure the best care for your child.
 - Text: I, [Parent/Guardian Name], consent to my child, [Child's Name], using the CleftCare app, which supports care management for cleft conditions by providing personalized resources, appointment tracking, and educational materials. I understand the app will collect and store necessary personal and health-related information, which will be securely handled in compliance with privacy regulations. I acknowledge that I can revoke this consent at any time by contacting CleftCare support at [support contact information].
 - Add Signature button
- Screen 3: Add Signature**
 - Header: < Add Signature
 - Instruction: Please provide additional details and then add your full name as signature.
 - Form field: Eli Powers (with person icon)
 - Signature area: A large square box containing a handwritten signature.
 - Buttons: Clear, Add Your Sign Here, Save
 - Add Signature button

All three screens share a common bottom navigation bar with icons for Home, List, and Profile.

Recording a Patient: On the recording screen, a central microphone button is used to start and stop each recording, while text displayed above the button presents the speech stimulus prompts which the CHW will be reading it to the child and then the child must repeat it. When the recording is started, the microphone button turns green and shows a text below it that the recording has started. When recording for a given stimulus is complete, the CHW taps the “Next” button in the top-right corner to save that audio to a secure server accessible to AIISH SLPs for later analysis, and this process is repeated for each stimulus in the list. “Next” will be enabled to go to the next phrase only when the current phrase is recorded once. After all stimuli have been recorded, the app displays a confirmation message stating “All audio is saved” to reassure the CHW that every sample has been securely uploaded.

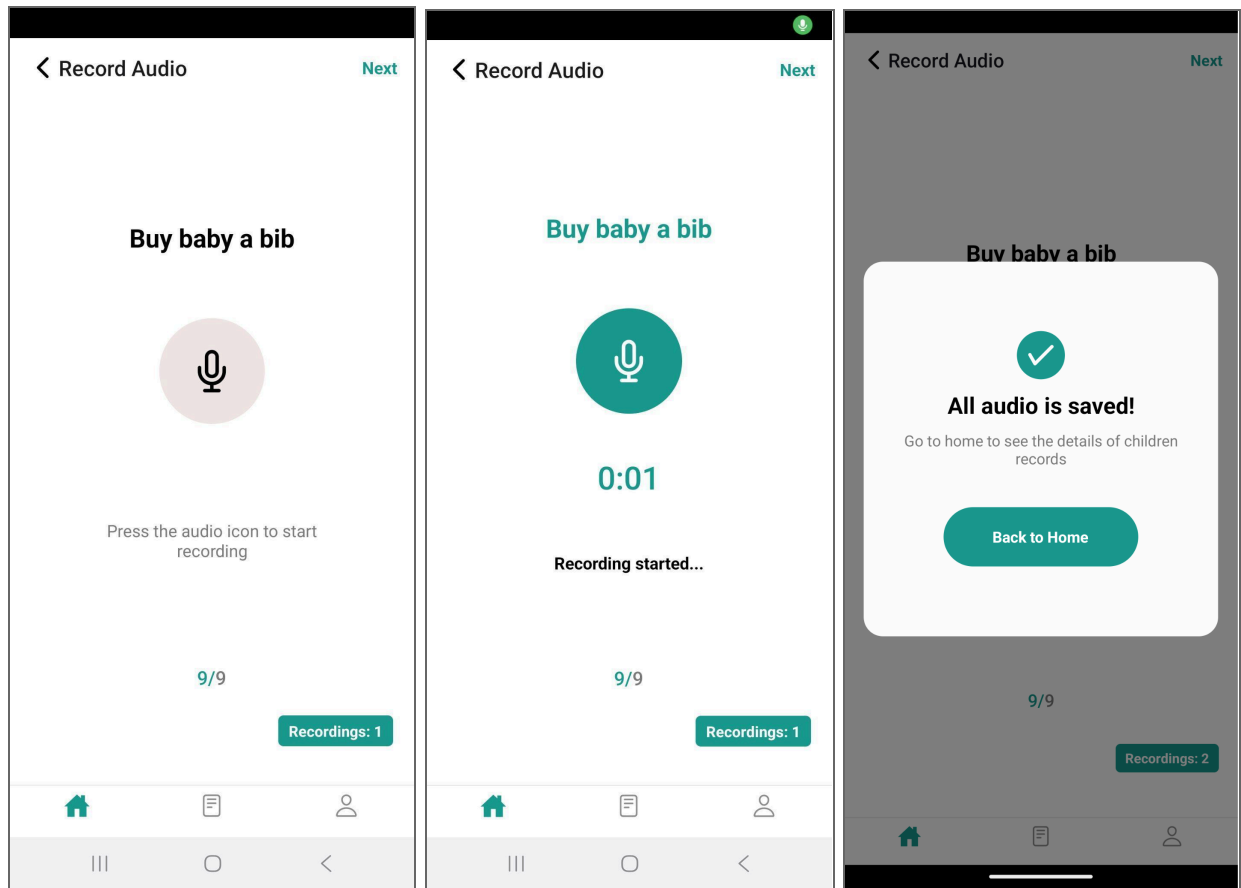

Edit Patient Details: The Home page lists all registered patients and allows CHWs to edit any patient's record. To edit a record, the CHW taps the desired patient entry to open their details and can update any demographic fields as well as their audio recordings, and then taps "Update" once all changes are complete. By selecting "Edit recordings," the CHW may re-record the phrases from start to end in one go and cannot skip to any other phrase before recording the current one at least once. The newly recorded audio samples will be automatically stored on the main server simultaneously.

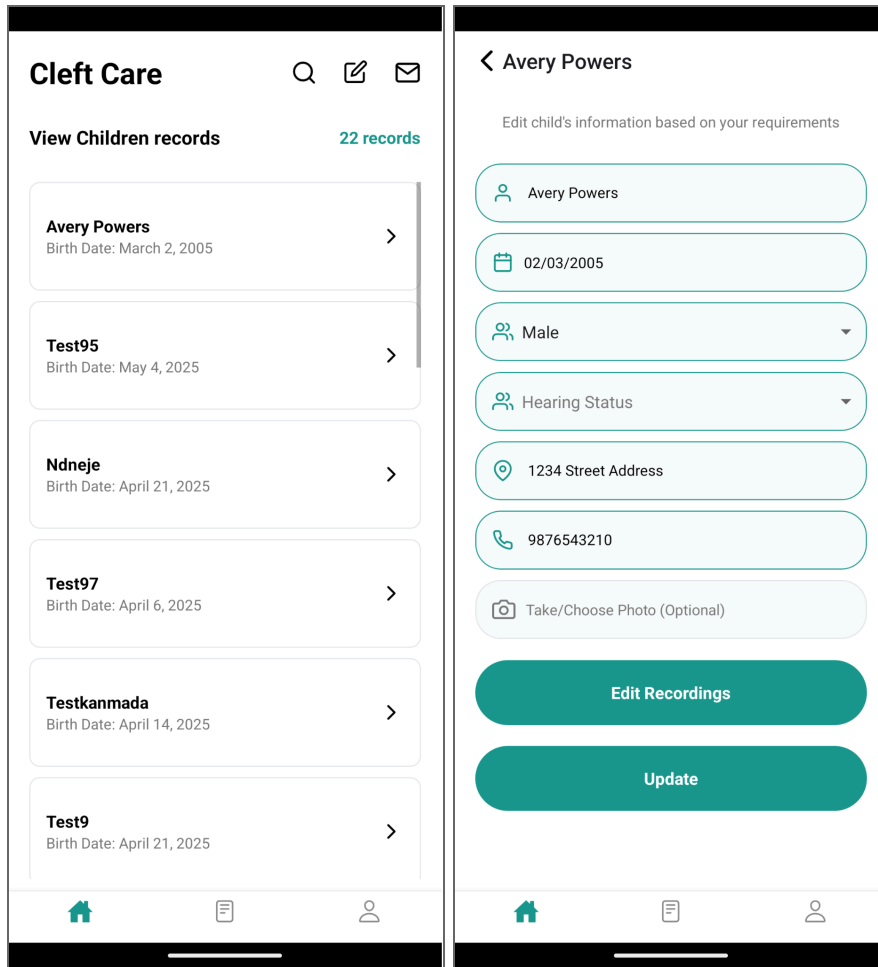

Help and Email Support: A Help button on the Home page (top-right) opens a simple support form where the CHW types their question or concern into an email-style message box and then taps “Submit” to send it to the clinician help center. The clinical team reviews these messages and responds to the CHW via email, addressing any technical or clinical queries related to the app or the recording process.

Help Center

×

Need help? Please let us know

cleftcareasu@gmail.com

Enter your message here

Submit

CHW Profile Page: The CHW Profile page displays the CHW's name, registered email address, and clinic location so they can quickly verify their account information. At the bottom of this page, a "Need Help?" text button takes the user directly to the Help page, where they can send an email to the clinician regarding any question or concern about using the app.

